EMPORIO
UNIQUE

EXCLUSIVE LINE

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Reflection therapy foot massage device E-3008

Manual

Thanks for choosing Reflection therapy foot massage device E-3008. Please read this manual carefully before using, especially for related attentions. Keep this manual properly for any necessary time.

Contents

Features of Product

Spare Parts and Function

Control Panel Instruction

Operation Instruction

Product Maintenance and Attentions


Simple Malfunction Analysis and How to Deal with

Specification

Features of Product

- Nice appearance, portable and convenient;
- Combination with shiatsu massage, kneading massage, air pressure massage and heating function together
- Kneading & shiatsu massage from toes to heels to cover all acupuncture points.
- Air pressure massage function makes you more relaxed.
- Six level air pressure massage intensities adjusted.
- Auto working mode and manual working mode available, two auto working modes adjusted.
- Two level heating temperatures adjustment: low, high temperature. 15 minutes auto shut off function.
- Power by DC12V adapter, safe and reliable.

Spare Parts and Function

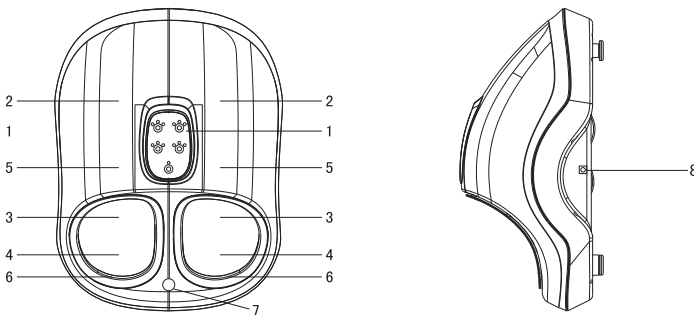1. Control panel: control the products' massage functions.
2. The position of rolling massage heads: massage forefoot area.
3. The position of shiatsu massage heads: massage the foot arch area.
4. The position of the kneading massage head: massage the heel area.
5. The position of the air bag: provide air pressure massage function.
6. Detachable and washable cloth cover.
7. DEFLATE BUTTON - use this button if you feel preassure too strong, hold this button for 5 seconds to release the air from bag.
8. DC12V/3A adapter.

Control Panel Instruction

1. **Kneading button:** selection of two levels massage speed
2. **Air Massage button:** Six level air pressure massage intensities selected in the manual working mode
3. **Auto button:** Two kinds of automatic working modes selected
4. **Heating button:** Low or high temperature selected in the manual or automatic working mode
5. **Power button:** turn on /off the massager

Operation Instruction

Plug output jack of adapter into power cord jack of the foot massager, then plug in-put jack of adapter into AC power socket, the buzzer makes a sound of “di”, the LED indicators of power flash, the foot massager is ready for working now.

Manual Working Mode

Manual working mode works only after the automatic working mode turns off, the below functions will work when manual working mode work.

1. Press “kneading” button - select two level massage speed slow/fast
2. Press “air massage” button button can select six level air massage intensities. (The first indicator light (marked with low) in green, it means the 1 st level. The first indicator light in red, it means the 2nd level. The second indicator light (marked with middle) in green, it means the 3rd level. The second indicator light in red, it means the 4th level. The third indicator light (marked with high) in green, it means the 5th level. The third indicator light in red, it means the 6th level.)
3. Press “Heating” button can select low, high temperature. (It can be adjusted in automatic working mode or manual working mode)

Automatic Working Mode

1. Turn on the power button, automatic mode A start to work(the indicator light on the left is on), the indicator lights of kneading function and air massage function are on, the below functions will work:

- A. Combination air pressure massage with kneading function
 - B. Air pressure massage function
 - C. Kneading massage function The above three kinds of function in cycle 3 times to the auto shut off time of 15 minutes to end
2. Press "Auto" button, automatic mode B start to work (the indicator light on the right is on), the below functions will work:
- A. Air pressure massage mode
 - B. Kneading massage mode
 - C. Combination air massage mode with kneading massage mode The above three kinds of function in cycle 3 times to the auto shut off time of 15 minutes to end
3. Press "Auto" button again, automatic working mode is off, and the indicator light is off.

This massager has the function of "15 minutes auto shut off", when the massage time is reach to 15mins, it stops working.

Product Maintenance and Attentions

For keep the massager in best condition, please follow up the correct usage and pay attention to below points.

1. Please do not massage the area on the foot where there is trauma or ulceration.
2. To prevent overload of the massage, please don't stand on it or sit on it.
3. At the beginning, please don't use the massage continuously for a long time. You could increase the using time after you adapting to it.
4. Please use it for a duration less than 30 minutes each time.

The person who is not sensitive of the heat needs to pay more attention to the surface of the appliance.

5. During the use, if any abnormal conditions spotted, such as overheating or smell, please stop using immediately and power off the massager.
6. After use, please power off the system by press "Power" button on the controller firstly, then disconnect the power plug from AC power. Please don't drag the power cord out forcibly.
7. Please unplug the adapter from the AC power jacket before cleaning the products.
8. Please use soft cotton sheet and neutral detergent to clean the massager
9. If you are receiving medical treatment or any of the following cases, please consult your doctor and use it in accordance with doctor's orders.
 - Cancer
 - Severe heart disease
 - Cardiac pacemaker
 - Pregnant women or children
10. Please don't drag or pull the power cord of the massager forcibly.
11. Please don't use the massage in water or any humid sites. If the massage falls in water accidentally, please turn off the power and unplug the power plug from socket immediately before fetching it out from water.
12. Please don't put the product under the pillow or the blanket when you using it.
13. Please do not use something to cover the massager when using it, if covering the massager, the motor temperature will be too high to cause malfunction or damage.
14. Please don't use the massager near stoves or any other heat source
15. Do use original accessories or part from the manufacturer of this massager
16. Please don't detach or dismantle the massage to repair by yourself without authorization from us. Any such behaviour will make the user unqualified to receive our guarantee along with after-sales service.
17. For your personal safety and life of this massager, please don't open the zipper. We will be not responsible for any consequences caused by such behaviour.

Malfunction Analysis and How to Deal With

Phenomena	Possible Cause	How to Deal With
Does not work?	1. The power supply is not on	1. Connect power supply
	2. Power plug is poor contact	2. Check the plug
	3. "Power" button has not been pressed	3. Press "power" button
Stop working?	1. Power off	1. Check whether the power supply interface is connected well
	2. Massage time is over 15 mins	2. Turn on "power" button again

Specification

Product name: Multifunctional Foot Massager
 Model No.: E-3008
 Voltage: AC100V~240V/50HZ~60HZ switch to DC12V/3.0A
 Power: 36W

1. This appliance does not intend to provide to the people who are physical disability, sensory or nerve defect or lack of experience and common sense to use, unless they got the supervision and guidance about how to use this equipment from the people who take responsibility to their safety.
2. Children must be supervised to ensure that they will not take instruments to play.